

Technical Data Sheet


PI/RNase Staining Buffer

Product Information

Material Number:	550825
Size:	100 ml
Storage Buffer:	Aqueous buffered solution containing $\leq 0.09\%$ sodium azide.

Description

Propidium Iodide (PI) is a fluorescent vital dye that stains DNA and RNA. PI binds to both DNA and RNA, so the latter must be removed by digestion with ribonucleases. The content of DNA as determined by flow cytometry, can reveal useful information about the cell cycle and the proteins involved in the regulation of the cell cycle. Cells in G2 and M phases of the cell cycle have double the DNA content of those in G0 and G1 phases. Cells in S phase have a DNA content lying between these extremes. PI is detected in the orange range of the spectrum using a 562-588 nm band pass filter. This reagent may be used to analyze cell cycle by flow cytometry in addition to use with antibodies for examining the expression of proteins during the cell cycle.

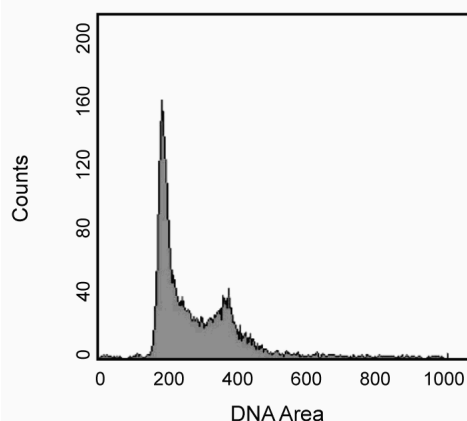

DNA content histogram of fixed Jurkat cells (Human T-cell leukemia; ATCC TIB-152). Jurkat cells were fixed with 1% paraformaldehyde (methanol free) and stored in 70% ethanol at -20°C . Cells were stained with 0.5 mL PI/RNase staining buffer for 15 minutes at room temperature and analyzed by flow cytometry.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.
Do not dilute or boil.

Application Notes

Application

Flow cytometry	Routinely Tested
----------------	------------------

Recommended Assay Procedure:

Flow cytometry: After fixing and permeabilizing your cell sample, use 0.5 mL /test (1×10^6 cells) and incubate for 15 minutes at room temperature before analysis. Please refer to http://www.bdbiosciences.com/documents/BD_FlowCytometry_DNA_Staining_Protocol.pdf for more protocol information.

Suggested Companion Products

Catalog Number	Name	Size	Clone
554656	Stain Buffer (FBS)	500 ml	(none)
554655	Fixation Buffer	100 ml	(none)
558052	Perm Buffer II	125 ml	(none)
558050	Perm Buffer III	125 ml	(none)

Product Notices

- Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- Avoid contact with skin and eyes.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



References

Douglas RS, Tarshis AD, Pletcher CH Jr, Nowell PC, Moore JS. A simplified method for the coordinate examination of apoptosis and surface phenotype of murine lymphocytes. *J Immunol Methods*. 1995; 188(2):219-228. (Biology)

Kalejta RF, Shenk T, Beavis AJ. Use of a membrane-localized green fluorescent protein allows simultaneous identification of transfected cells and cell cycle analysis by flow cytometry. *Cytometry*. 1997; 29(4):286-291. (Biology)

BD Biosciences

bdbiosciences.com					
United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be constructed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD

